Heparan Sulfate Glycan Array User Manual



Website: <u>http://www.zbiotech.com/home.html</u> Tel: (720) 285-3587 Email: <u>info@zbiotech.com</u>

For Research Use Only Copyright 2021, Z Biotech, LLC. All Rights Reserved.

Introduction

Heparan sulfate is a sulfated, unbranched polysaccharide found on cell surfaces as well as in the extracellular matrix. Heparan sulfates interact with a variety of proteins and are known to be important to the regulation of many physiological and pathophysiological functions. Because of the prevalence and multifunctionality of heparan sulfates, investigation of the binding roles of this type of glycan has become invaluable for understanding how to combat related diseases.

Z Biotech has expanded its microarray technology to provide a research platform specifically for the investigation of heparan sulfates. Our Heparan Sulfate Array can be used to determine the specific binding interactions between various types of heparan sulfates and their potential protein partners via microarray assays.

This manual is provided as a comprehensive guide to help the researcher acquire clear results from the assay. Please read through carefully before starting your experiment.


Handling and Storage

Store the bag of slides and any buffers in a 4°C refrigerator if they are to be assayed within 24 hours upon receipt. For long term storage keep the bag of slides at -20°C. Avoid freezing and thawing multiple times. Purchased slides and buffers should be used within 6 months.

Allow the bag of slides to equilibrate to room temperature at least 20 minutes before opening. After opening, re-seal any unused slides in the moisture barrier bag with a desiccant inside and refreeze.

Array Map/Schematic

24 Heparan Sulfate Glycan Array slides have either 8 or 16 subarrays. Arrays are printed on the side with the "Z Biotech" label and 4-digit number ID facing upward. The "Z Biotech" label is located on the bottom center from a landscape view. The number ID is consistent with the barcode ID on the bottom from a portrait view. Dimensions and array maps are shown below.

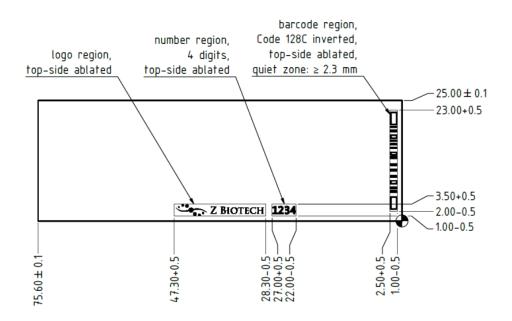


Array Map (16-subarray slides):

ide	
ray Sl	Sub- array Sub- sub- array Sub- array sub- array Sub- array Sub- array
-subarray Slid	ub- uub- trray trub- trub- trub- trub- trub- trub- trub- trub- trub-
16-s	

HS001	HS001	HS001	HS001	HS001	HS001	HS002	HS002	HS002	HS002	HS002	HS002
HS003	HS003	HS003	HS003	HS003	HS003	HS004	HS004	HS004	HS004	HS004	HS004
HS005	HS005	HS005	HS005	HS005	HS005	HS006	HS006	HS006	HS006	HS006	HS006
HS007	HS007	HS007	HS007	HS007	HS007	HS008	HS008	HS008	HS008	HS008	HS008
HS009	HS009	HS009	HS009	HS009	HS009	HS010	HS010	HS010	HS010	HS010	HS010
HS011	HS011	HS011	HS011	HS011	HS011	HS012	HS012	HS012	HS012	HS012	HS012
HS013	HS013	HS013	HS013	HS013	HS013	HS014	HS014	HS014	HS014	HS014	HS014
HS015	HS015	HS015	HS015	HS015	HS015	HS016	HS016	HS016	HS016	HS016	HS016
HS017	HS017	HS017	HS017	HS017	HS017	HS018	HS018	HS018	HS018	HS018	HS018
HS019	HS019	HS019	HS019	HS019	HS019	HS020	HS020	HS020	HS020	HS020	HS020
HS021	HS021	HS021	HS021	HS021	HS021	HS022	HS022	HS022	HS022	HS022	HS022
HS023	HS023	HS023	HS023	HS023	HS023	HS024	HS024	HS024	HS024	HS024	HS024
NC	NC	NC	NC	NC	NC	PC1	PC1	PC1	PC1	PC1	PC1
PC2	PC2	PC2	PC2	PC2	PC2	PC3	PC3	PC3	PC3	PC3	PC3
PC4	PC4	PC4	PC4	PC4	PC4	Marker	Marker	Marker	Marker	Marker	Marker

Array Map (8-subarray slides):

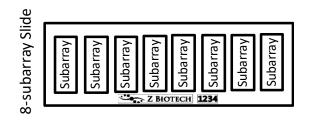

HS001	HS001	HS001	HS001	HS001	HS001	HS002	HS002	HS002	HS002	HS002	HS002	HS003	HS003	HS003	HS003	HS003	HS003
HS004	HS004	HS004	HS004	HS004	HS004	HS005	HS005	HS005	HS005	HS005	HS005	HS006	HS006	HS006	HS006	HS006	HS006
HS007	HS007	HS007	HS007	HS007	HS007	HS008	HS008	HS008	HS008	HS008	HS008	HS009	HS009	HS009	HS009	HS009	HS009
HS010	HS010	HS010	HS010	HS010	HS010	HS011	HS011	HS011	HS011	HS011	HS011	HS012	HS012	HS012	HS012	HS012	HS012
HS013	HS013	HS013	HS013	HS013	HS013	HS014	HS014	HS014	HS014	HS014	HS014	HS015	HS015	HS015	HS015	HS015	HS015
HS016	HS016	HS016	HS016	HS016	HS016	HS017	HS017	HS017	HS017	HS017	HS017	HS018	HS018	HS018	HS018	HS018	HS018
HS019	HS019	HS019	HS019	HS019	HS019	HS020	HS020	HS020	HS020	HS020	HS020	HS021	HS021	HS021	HS021	HS021	HS021
HS022	HS022	HS022	HS022	HS022	HS022	HS023	HS023	HS023	HS023	HS023	HS023	HS024	HS024	HS024	HS024	HS024	HS024
NC	NC	NC	NC	NC	NC	PC1	PC1	PC1	PC1	PC1	PC1	PC2	PC2	PC2	PC2	PC2	PC2
PC3	PC3	PC3	PC3	PC3	PC3	PC4	PC4	PC4	PC4	PC4	PC4	Marker	Marker	Marker	Marker	Marker	Marker

No.	Structure
HS001	GlcNAc-GlcA-GlcNAc-GlcA
HS002	GlcA-GlcNAc-GlcA-GlcAA
HS003	GIcNAc-GIcA-GIcNAc-GIcA-GIcNAc-GIcA
HS004	GICA-GICNAc-GICA-GICNAc-GICAAGICNAC-GICA
HS005	GlcNAc-GlcA-GlcNAc-GlcA-GlcNAc-GlcA-GlcNAc-GlcA
HS006	GlcA-GlcNAc-GlcA-GlcNAc-GlcA-GlcNAc-GlcA-GlcNAc-GlcA
HS007	GlcNS-GlcA-GlcNS-GlcA
HS008	GlcA-GlcNS-GlcA-GlcNS-GlcA
HS009	GlcNS-GlcA-GlcNS-GlcA-GlcNS-GlcA
HS010	GIcA-GIcNS-GIcA-GIcNS-GIcA-GIcNS-GIcA
HS011	GlcNS-GlcA-GlcNS-GlcA-GlcNS-GlcA-GlcNS-GlcA
HS012	GIcA-GIcNS-GIcA-GIcNS-GIcA-GIcNS-GIcA-GIcNS-GIcA
HS013	GIcA-GIcNS-GIcA-GIcNS-GIcA-GIcNS-GIcA-GIcNS6S-GIcA
HS014	GlcA-GlcNS-GlcA-GlcNS-GlcA-GlcNS6S-GlcA-GlcNS6S-GlcA
HS015	GlcA-GlcNS-GlcA-GlcNS6S-GlcA-GlcNS6S-GlcA-GlcNS6S-GlcA
HS016	GlcA-GlcNS6S-GlcA-GlcNS6S-GlcA-GlcNS6S-GlcA-GlcNS6S-GlcA
HS017	GlcNS-GlcA-GlcNS-GlcA-GlcNS-IdoA2S-GlcNS-GlcA
HS018	GlcNS-GlcA-GlcNS-IdoA2S-GlcNS-IdoA2S-GlcNS-GlcA
HS019	GlcNAc-GlcA-GlcNS-IdoA2S-GlcNS-IdoA2S-GlcNS-GlcA
HS020	GlcNS6S-GlcA-GlcNS6S-GlcA-GlcNS6S-IdoA2S-GlcNS6S-GlcA
HS021	GlcNS6S-GlcA-GlcNS6S-IdoA2S-GlcNS6S-IdoA2S-GlcNS6S-GlcA
HS022	GlcNAc6S-GlcA-GlcNS6S-IdoA2S-GlcNS6S-IdoA2S-GlcNS6S-GlcA
HS023	GlcNS6S-GlcA-GlcNS3S6S-IdoA2S-GlcNS6S-GlcA
HS024	GlcNAc6S-GlcA-GlcNS3S6S-IdoA2S-GlcNS6S-IdoA2S-GlcNS6S-GlcA

24 Heparan Sulfate Glycan Identification List:

Controls

- NC: Negative control, Print Buffer
- PC1: Positive control 1, Biotinylated Mannose (0.01 mg/mL)
- PC2: Positive control 2, Human IgG (0.1 mg/mL)
- PC3: Positive control 3, Mouse IgG (0.1 mg/mL)
- PC4: Positive control 4, Rabbit IgG (0.1 mg/mL)

Marker: Anti-human IgG, Cy3 (0.01 mg/mL) and anti-Human IgG, Alexa647 (0.01 mg/mL)

Materials Required

- Arrayed glass slides
- 16 or 8 cassettes
- Glycan Array Blocking Buffer (GABB, Item #10106), add 1% BSA (10 mg/ml) if needed
- Glycan Array Assay Buffer (GAAB, Item #10107), add 1% BSA (10 mg/ml) if needed
- Wash Buffer: 50 mM Tris-HCl, 137 mM NaCl, 0.05% Tween 20, pH 7.6
- Laser fluorescence scanner (able to scan at the wavelength of your fluorophore)
- Coplin jar
- Adhesive slide cover film

Preparation of assay samples:

Prepare glycan-binding protein samples or detection antibodies in a centrifuge tube by diluting with the GAAB. For the fluorescently labelled streptavidin we recommend a concentration of 1 μ g/mL. We suggest a range of 50 μ g/ml to 0.1 μ g/ml concentration for protein samples, although some experimentation may be required to establish the concentration that will provide the highest binding signals with the lowest background fluorescence. This can be accomplished by applying a different dilution of samples to different wells of the array. In addition to testing a dilution range for your protein of interest, we recommend setting up control assays for any additional detection or secondary antibodies to ensure that any binding observed is specific to your protein of interest. This may be especially important as heparin sulfates have been known to show weak non-specific binding to proteins. A fluorescent signal due to specific binding to your protein of interest should be dose-dependent within the dynamic range of your protein dilution, and should have a positive binding signal after a signal from control assays has been subtracted. Calculate the volume of sample needed depending on how many slides and subarrays are to be assayed. Use 100 µL volume of sample per well for 16 subarray cassettes and 200 µL for 8 subarray cassettes to ensure full and even coverage of the printed area throughout incubation. If necessary, the assay can be done successfully with a minimal volume of 60 μ L per well for 16 subarray cassettes and 80 μ L for 8 subarrays. Using a minimal volume in the wells has an increased risk of the array drying out during the assay, and may also cause unequal distribution of the sample across the arrayed surface which may result in signal variation. Please ensure each sample is homogeneous and thoroughly mixed.

Assay Protocol

Part 1 - Blocking

Handle the slide in a clean, dry environment. Use gloves and avoid touching the slide surface

- 1. Allow the arrayed slides to equilibrate to room temperature (20-30 minutes) before opening the moisture barrier bag.
- 2. Add blocking buffer to each subarray well. We recommend using 100 µL per well for 16 subarray cassettes and 200 µL for 8 subarray cassettes to ensure full and even coverage of the printed area throughout incubation.
- 3. Cover the wells with adhesive film to prevent evaporation and incubate slide on shaker at 85 rpm for 1 hour. Longer incubation time is acceptable, but not necessary.

Make sure the orbital shaker is completely flat. If the slide is sloped in any direction during incubation it can cause variation in binding.

Part 2 – Binding assay

1. Unless the glycan-binding sample of interest is bacteria or cells, centrifuge samples briefly to avoid adding irrelevant particles to the array.

- 2. Touch the pipette tip to the corner of the well of the cassette and tip the slide so that the sample pools to that corner to remove blocking buffer. Avoid touching the array surface.
- 3. Immediately apply the glycan-binding sample of interest to each well. We recommend using 100 μL per well for 16 subarray cassettes and 200 μL for 8 subarray cassettes to ensure full and even coverage of the printed area throughout incubation. Avoid leaving air bubbles.
- 4. Seal the wells with adhesive film to prevent evaporation. If the sample is fluorescently labelled, cover with aluminum foil to keep it in the dark. Incubate on the shaker for 1-3 hours at 100 rpm. If the samples can easily aggregate, shake at higher speed to prevent aggregation. Longer incubation time may increase binding signal, especially for weakly binding samples.

Avoid allowing the slides to dry out at any point during the assay, especially during long incubation times. Make sure the adhesive film is sealed around each well.

If your glycan-binding samples are fluorescently labelled, go directly to Part 6 – Final wash and dry.

Part 3 - Wash

- 1. Remove glycan-binding samples from each well by gently touching the pipette tip to the corner of the well of the cassette and tip the slide so that the sample pools to that corner. Avoid touching the array surface, but a gentle touch is okay to ensure no sample is left pooled in the corners.
- 2. Add wash buffer to each well. We recommend using 100 µL per well for 16 subarray cassettes and 200 µL for 8 subarray cassettes. Cover the wells with adhesive film and incubate on the shaker for 5 minutes at 85 rpm. Completely remove the wash buffer by pipette and repeat this step. Avoid allowing the slide to dry out and have your next wash or sample ready before you remove the wash buffer.

If your glycan-binding sample is biotinylated, go directly to Part 5 – Fluorescent staining.

Part 4 - Binding of biotinylated antibody (Sandwich Assay Format)

- 1. Unless the secondary biotinylated antibody sample is bacteria or cells, centrifuge samples briefly to avoid adding irrelevant particles to the array.
- 2. After completely removing the wash buffer immediately add the biotinylated antibody to each well. We recommend using 100 μ L per well for 16 subarray cassettes and 200 μ L for 8 subarray cassettes. Seal the wells with adhesive film and incubate on the shaker for 1 hour at 100 rpm. Longer incubation time is acceptable, but not necessary.
- 3. After incubation repeat Part 3 Wash

Part 5 - Fluorescent staining

- 1. Centrifuge fluorescent-labeled streptavidin samples briefly to avoid adding irrelevant particles to the array.
- 2. After completely removing the wash buffer immediately add the fluorescently labelled streptavidin sample. 100 μ L per well is recommended for 16 subarray cassettes and 200 μ L for 8 subarray cassettes. Seal the wells with adhesive film and shield the wells from light with aluminum foil. Incubate on the shaker at 85 rpm for 1 hour. Longer incubation time is acceptable, but not necessary.

Part 6 – Final wash and dry

1. Touch the pipette tip to the corner of the well of the cassette and tip the slide so that the sample pools to that corner to remove it. Avoid touching the array surface.

- Briefly rinse each well with wash buffer. 100 μL per well is recommended for 16 subarray cassettes and 200 μL for 8 subarray cassettes.
- 3. Completely remove the wash buffer by pipette. Avoid touching the array surface. Repeat steps 2 and 3.
- 4. Disassemble the cassette from the slide. For the provided cassette this can be done by holding the slide with one hand at the top and bottom edges and sliding out the cassette clips one by one with the other hand. If your provided cassette has metal clips, they can be removed by rotating the clip outwards from the bottom of the slide. When the clips have been removed place the slide on the table and hold a small outer edge of the slide to the table as you gently peel the cassette off.
- 5. Immediately immerse the slide in a coplin jar or beaker full of wash buffer. Do not touch the surface of the array or allow the array surface to touch the sides of the beaker or jar.
- 6. Place the jar or beaker on the 60 rpm shaker for 10 minutes.
- 7. Decant the wash buffer from the jar or beaker while holding the slide in place (only touch the edge of the slide) and then add sterile de-ionized water to immerse the slide.
- 8. Place the jar or beaker on the 60 rpm shaker for 2 minutes.
- 9. Decant the water from the jar or beaker.
- 10. Allow the slide to dry completely in a clean, dust free environment before scanning.

Analysis

Scan the slide in a laser fluorescence scanner at the wavelength of emission for the fluorophore used. Adjust the laser power and PMT to obtain the highest possible signals without any being saturated. Analyze data with microarray analysis software. If there is specific binding the signal intensity should be significantly higher than the background signal (area where there are no printed spots). Fluorescent signal due to specific binding to your protein of interest should be both dose-dependent with your protein dilution (unless the protein concentration range is too high and saturates the spots), and should have positive binding signal after signal from control assays has been subtracted. Our standard method of comparing signal intensities is to quantify the median signal intensity data and subtract the background intensity. Subtracting signal from negative control spots as well as the same spots on a negative control assay will give more accurate specific binding data.

Interpretation of Control Signals:

<u>Negative Control (Print Buffer)</u>: The negative control should produce a signal close to the intensity of the background. Since there is no binding involved with the negative control, any other signals around the negative control's intensity are also not binding. Subtracting the negative control's signal from the other binding signals will give more accurate specific binding data.


<u>Biotinylated Mannose (PC1)</u>: This positive control will bind directly to the fluorescent labelled streptavidin. If your glycanbinding protein sample is already fluorescently labelled, or in any case where the addition of fluorescent labelled streptavidin to the array was not preformed (Part 5 – Fluorescent staining) this positive control will not be reactive.

IgG (PC2, PC3, PC4): IgG is an antibody found in blood that is a primary component of humoral immunity. If the glycanbinding or secondary antibody sample is an anti-IgG from human, rabbit, or mouse it should bind to the respective IgG control.

<u>Marker</u>: The array marker should show a fluorescence signal regardless of the assay. It is there primarily to aid with orientation of the array map during analysis.

Typical Binding Assay Result from the 24 Heparan Sulfate Glycan Array

Example: 24 Heparan sulfate glycan array on 16 subarray formats. A subarray assayed with protein antithrombin III (10 μ g/ml), followed by anti-antithrombin III antibody (2 μ g/ml), and then anti-mouse IgG antibody AlexaFluor555 (10 μ g/ml). The array was scanned with GenePix scanner at 500 PMT and 100% laser power at 532nm wavelength. Positive controls 2-4 and the marker show binding as expected with the mouse IgG control being the strongest (PC3). HS023 and HS024 show specific binding to Antithrombin III.

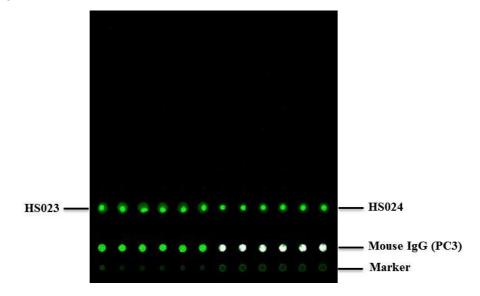

Troubleshooting

Condition	Possible Causes	Potential Solutions					
High Background	 Concentration of glycan-binding protein samples is too high Concentration of fluorescent samples is too high. Arrays are not thoroughly washed Slide drying out during assay Excessive particles in the samples due to sample aggregation, dust, etc. 	 Use a lower concentration range of samples. Consider a wider range if you are unsure where the detection limit is. Apply longer times for washing steps and use a higher shaking rate Make sure wash buffer and sample is completely removed before the next step Make sure adhesive film fully seals the wells to avoid evaporation Centrifuge the samples prior to assay to avoid adding irrelevant particles If you think that the protein is aggregating during incubation, try shaking at a higher speed 					
Signal Variation	 Slide drying out during assay. Binding samples are not equally distributed in the wells Glycan-binding protein aggregation during incubation Bubbles during incubation 	 Make sure wells are sealed to prevent evaporation during incubation. Apply a larger volume of sample to each well to ensure equal distribution Use a higher shaking rate during incubation Make sure samples are homogeneous, mixed thoroughly, and do not leave bubbles on the array surface 					
Unexpected Binding	 Cross contamination between wells or other sources. Sample contamination 	 Make sure to use sterilized pipette tips and tubes used for sample application and preparation Ensure cassette is pressed firmly to the slide so that there are no gaps to allow leaking between wells Be careful not to cross contaminate samples when applying to the wells, even during wash steps 					